



Technical specifications

Bed height

Low position – surface edge to floor.....18" (45.7 cm)
 with optional IntelliDrive® power transport.....22" (55.9 cm)
 Maintains low height when placed in chair.....Yes
 High position – surface center to floor.....37" (94.0 cm)
 Headboard/push handles maintain same height regardless of patient
 surface height.....Yes

Sleep surface

Surface thickness.....8" (20.3 cm)
 Width.....35.5" (90.2 cm)
 Length, fully extended.....86" (218.0 cm)
 Length, fully retracted.....75" (191.0 cm)

Overall width

Siderails stowed.....37" (94.0 cm)
 Siderails up.....40" (102.0 cm)

Overall length

Fully extended.....94.5" (240.0 cm)
 Fully retracted.....82.5" (210.0 cm)

Bed angles

Head section.....0°-65°
 Knee gatch.....0°-16°
 Foot section.....0°-27°
 Max Trend/reverse Trend.....15°/10°
 Head and Trend angle indicators.....Yes

Caster & braking system

Optional IntelliDrive® power transport.....Yes
 Caster diameter.....5" (12.7 cm)
 with IntelliDrive® power transport.....6" (15.2 cm)
 Four-wheel braking.....Yes
 Brakes lock both rolling and swiveling.....Yes

Available surfaces

Comfortline® prevention surface.....Standard
 PrimeAire® ARS pressure relief surface.....Optional
 Active Integrated Response™ (A.I.R.™) treatment surface.....Optional

A.I.R. surface standard features:

- Turn Assist
- Three active, independently-controlled zones
- Max Inflate
- Automatic pressure adjustment for weight and bed position

Features common to all surfaces:

FlexAfoot™ retractable foot mechanism.....Standard
 Shear liner.....Standard
 Firm, grippable edge.....Standard
 Anti-hamocking surface.....Standard
 Meets all applicable flammability standards.....Yes

Weight limit (safe working load).....500 lbs. (227 kg)

Optional bed exit alarm system

Sensitivity levels.....3
 Modes.....Positioning, Exiting, Out-of-Bed

Optional scale system

Accuracy.....1% of patient weight
 Weigh in any position (except Trend/Rev Trend).....Yes
 Scale weight capacity.....500 lbs. (227 kg)
 (OIML compliant scales available in Europe)

Patient controls

Siderail patient bed controls.....Standard
 Docking pendant.....Optional
 SideCom® patient communication controls
 (Nurse call, universal TV/radio, and lighting).....Optional

Electrical characteristics

Backup battery power.....Standard
 Safety listings.....UL2601 & IBC 60601-1



At a glance

Improve patient safety and satisfaction:

- Lowest height of any acute care bed helps reduce the risk and impact of patient falls.
- Firm, supported perimeter of surface enables safer patient ingress and egress.
- Minimized siderail gaps reduce opportunities for entrapment.
- Help patients achieve greater comfort in LowChair® position.

Protect your caregiving staff:

- Create a firm surface for easier patient transfers with Max Inflate.
- Perform patient examinations and linen changes more easily with the Turn Assist function in the VersaCare A.I.R.™ surface.
- Transport the patient and the bed more easily with the optional IntelliDrive® power-assisted transport feature.

Continuously deliver pressure relief:

- Provide immediate intervention with the "always on" wound therapy.
- Surface thicknesses ranging from seven to eight inches surpass other surface thicknesses for patient comfort.

The bed system without compromise

In a Med-Surg environment, many nurses want their patient care beds to help them deliver effective wound therapy, improve patient safety, and improve caregiver safety. In the past, it was difficult to find a comprehensive solution to these needs in one integrated system, but the wait is over. The new flexible and versatile VersaCare® bed system from Hill-Rom eliminates the need to compromise. It is designed to help achieve all three goals in a single bed choice, with the quality and comprehensive support expected from Hill-Rom.

The VersaCare bed is a full-featured patient care bed system that offers the lowest height in acute care, an integrated air surface that helps automate the delivery of wound care therapy, and many features to reduce the risk of caregiver injury. In addition, this latest innovation from Hill-Rom offers superior patient comfort and intuitive controls to help caregivers effectively and easily operate the system to its fullest capacity.



The new standard for patient safety

Med-Surg units include patients with a wide variety of needs. Some are awake and mobile, sitting up and getting out of bed with little or no assistance. Others need more help. The VersaCare bed is designed to help prevent or reduce the impact of falls for patients at a range of acuity levels.

Low height

For more active patients, the VersaCare bed's lower position—the **lowest height of any acute care bed**—means it's safer to get in and out of bed. For higher acuity patients, the VersaCare bed makes it easy to position patients up in chair without leaving the bed. And with the **LowChair®** feature, caregivers can keep the bed at the same low height while placing patients in a chair position, without compromising on the patient's freedom of movement or the caregiver's time.

Patient safety features

The VersaCare bed's unique surface shape supports the patient as he or she sits on the edge of the surface, providing firm support for egress. **Nightlights** at the side of the bed are focused on where the patient's feet will land. Comfortable, **easy-to-grip siderails** protect without restraining. Even the optional **Bed Exit Alarm System** has three levels of sensitivity, to help monitor without overly confining a patient's movement.



Bed plugs into hospital network to monitor unsafe bed conditions

The demands on caregivers increase every day, and Hill-Rom is leveraging computer network technology to enhance patient safety. VersaCare with OnSite™ technology can send data points from the bed to the hospital network for viewing, reporting, and alerting, using most existing computer hardware.



Maximum protection and efficiency for caregivers

Keeping caregivers safe and on the job is a concern of every care facility. The VersaCare bed protects caregivers by helping them avoid or eliminate risky tasks with these features:

LowChair® Position

The LowChair height enables caregivers to easily slide patient from bed to wheelchair. Even in the LowChair position, shorter patients can put their feet firmly on the floor.



Built-in Scale

This optional feature eliminates the need to move patients to a separate scale.

FlexAfoot™ Mechanism

Helps keep patients from sliding down in bed, reducing the need for caregivers to pull patients up.



Caregiver safety and ease

In addition to helping caregivers do their jobs safely, the VersaCare bed helps make routine tasks easier with these features:

IntelliDrive® Powered Transport System

An optional power-assist transport system, allows a single caregiver to easily transport patients throughout the facility.

Fixed headboard

The headboard stays at the optimal height for pushing, no matter how low the bed is. Patients can stay at a safer height without caregiver back strain.


Turn Assist

Available with the VersaCare A.I.R.™ surface, this feature makes it easier for a single caregiver to perform patient examinations, linen changes, and other routine tasks.



HandsFree® Foot Controls

Keep hands free to manage lines or guide a patient as the bed height or head angle is adjusted.

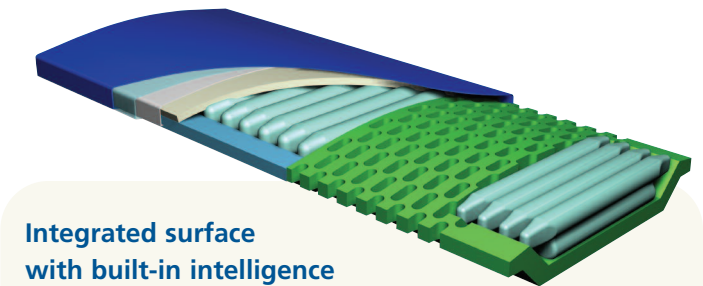

Integrated surface with built-in intelligence

Unlike other beds in the Med-Surg environment, the VersaCare® bed combines low height with therapy. Choose from our Active Integrated Response (A.I.R.™) treatment surface, our PrimeAire® ARS pressure relief surface, or our Comfortline® prevention foam surface. Because they are integrated into the bed design, they provide the best protection from falls and entrapment for patients. With an integrated A.I.R. and ARS surface, therapy is always being provided whenever the bed is plugged in. Unlike non-integrated therapy surfaces, there are no cumbersome, noisy boxes and no extra plugs, cords or hoses to trip over.

The VersaCare bed's surface choices take into account the main factors that contribute to decubitus ulcers:



Pressure

Working with the intelligence of the bed to accommodate patient weight and bed position, our A.I.R. treatment surface uses zones of continuous low pressure to reduce the peak pressures that cause and aggravate skin ulcers. PrimeAire ARS offers a non-powered air surface with low interface pressures to help achieve desirable clinical outcomes.



Temperature and moisture

Our A.I.R. treatment surface is constructed to assist in reducing the amount of sweat and moisture developed on the skin.



Shear

All of our surfaces have a shear liner and stretchable ticking to minimize shear on the patient's skin. With the Max Inflate function on our A.I.R. surface, transfers are easier for both patients and caregivers. FlexAfoot™ feature provides a customized surface to maintain proper patient positioning, thus helping to reduce shearing forces.