

TABLE OF CONTENTS

Refrigerators / Freezers	2
Accucold 220 lb Ice & Water Dispenser with Filter	3
1 Cu.Ft. Compact Medication/Vaccine Refrigerator, Certified to NSF/ANSI 456 Vaccine Storage Tabletop	3
Accucold Product USB Digital Data Logger (DDL)	4
ECOLAB OR SOLUTIONS ORS-1075HS/ORS-2038D HUSH SLUSH SYSTEM	5
Accucold 24" Wide All-Freezer up to -30 deg C capable	5
Accucold 3.08 Cu.Ft. -85°C Ultra Low Chest Freezer 36.25in H x 23.75in W x 28.0in D	6
Accucold 18 Cu.Ft. Upright Vaccine Refrigerator, Certified to NSF/ANSI 456 Vaccine Storage Standard	6
Accucold Under Counter 6 Cu.Ft. Vaccine Refrigerator, ADA Height	7
24" Wide All-Refrigerator/All-Freezer Combination	8
Accucold 160 lb Ice & Water Dispenser without Filter	9
Accucold 160 lb Ice & Water Dispenser with Filter	9
Accucold 450 lb Ice & Water Dispenser without Filter	10
Accucold 450 lb Ice & Water Dispenser with Filter	10
Accucold 4.8 Cu.Ft. -85°C Ultra Low Chest Freezer 36.25in H x 28.5in W x 28.0in D	11
Accucold 4.8 Cu.Ft. -85°C Ultra Low Chest Freezer 36.25in H x 28.5in W x 28.0in D with 9 Stainless Vertical Box Shelves	12
12 lb. Drain-Free Countertop Icemaker Cabinet Style	12
Summit 24" Wide All-Refrigerator, ADA Compliant	13
Accucold Ice & Water Dispenser with Filter	13
Accucold Under Counter Medication/Vaccine Refrigerator (Glass Door)	13
Accucold Under Counter Medication/Vaccine Refrigerator	14
Accucold 49 Cu.Ft. Upright Pharmacy Refrigerator	14
Accucold 24" Wide Stacked combination offers separate refrigerator/freezer storage	15
Accucold 49 Cu.Ft. Upright Pharmacy Freezer	16



REFRIGERATORS / FREEZERS

Auxo Medical offers a wide range of medical refrigerators and freezers designed to meet the specific needs of healthcare facilities. Our inventory includes both new and refurbished equipment from reputable brands like Accucold and Summit. These products are essential for storing vaccines, medications, and other temperature-sensitive materials, ensuring their efficacy and safety.

The refrigerators and freezers offered by Auxo Medical come in various sizes and configurations, from compact under-counter units, large upright models and ultra low (-30° c) Freezers. They feature advanced temperature control systems, ensuring precise and consistent temperature maintenance. Many models are also equipped with alarms that alert users to potential temperature fluctuations, preventing spoilage and ensuring product integrity.

The benefits of these refrigerators and freezers include reliable performance, energy efficiency, and compliance with industry standards. These products are designed to meet the demanding requirements of healthcare environments, providing peace of mind for healthcare professionals. Some of the specific uses for these products include vaccine storage, medication storage, blood storage, and laboratory sample storage.

Our catalog consists of brand new, reconditioned and refurbished equipment. Would you like to filter the equipment by:

[BRAND NEW](#)

[RECONDITIONED](#)

[REFURBISHED](#)



ACCUCOLD 220 LB ICE & WATER DISPENSER WITH FILTER

The **Accucold AIWD220 Ice & Water Dispenser** is a commercially approved countertop ice and water dispenser that fits slim spaces on standard counters. It's ETL-listed (NSF-12) with a stainless steel cabinet, black dispenser panel, and UV lighting in the dispenser for easy cleaning. Capable of up to 220 lbs of ice per day, it offers three dispensing options (ice, water, or ice+water) and adjustable legs for flexible installation. Top-vented with a drain requirement, it includes a removable tray and a drainage-ready design, making it ideal for hospital waiting areas, break rooms, and other high-use environments. Features filtration options are available separately to optimize performance. Specifications highlight a compact footprint (13.88" W x 23.75" D x 33.5" H), 115 V/60 Hz power, stainless steel front, and durable 1-year parts & labor warranty with a 5-year compressor warranty.

[Read More](#)

SKU: AIWD220FLTR



1 CU.FT. COMPACT MEDICATION/VACCINE REFRIGERATOR, CERTIFIED TO NSF/ANSI 456 VACCINE STORAGE TABLETOP

The **Accucold 1 Cu.Ft. Compact Vaccine Refrigerator** stands out as a reliable and efficient solution for vaccine storage, certified to the NSF/ANSI 456 Vaccine Storage Standard. Specifically designed for pharmacy, medication, and vaccine applications, this refrigerator supports adherence to CDC/VFC guidelines. It features a purpose-built design that ensures temperature stability with a $\pm 1^{\circ}\text{C}$ variation, monitored by an NTC sensor. The refrigerator includes a 2-year calibration certificate for record-keeping, while its antimicrobial silver ion-coated handle adds a layer of protection against germ spread. Users can adjust the interior temperature between $+2$ to $+8^{\circ}\text{C}$ using intuitive controls, and a suite of alarms—temperature, open door, power failure, and sensor failure—ensure optimal storage conditions. Additional features include remote alarm contacts, automatic defrost with adjustable cycle defrost, a factory-installed lock for security, and a microprocessor control panel with a resettable min/max display. With a compact design and a capacity of 1.0 cu.ft., the refrigerator is both practical and versatile, offering three adjustable wire shelves and an interior light for convenience. Its robust construction and advanced features make it an essential tool for ensuring the safe storage of vaccines and other sensitive medical materials.

[Read More](#)

SKU: ARG1PV456



ACCUCOLD PRODUCT USB DIGITAL DATA LOGGER (DDL)

The **Accucold USB Digital Data Logger (DDL)**, a cutting-edge solution for precise temperature monitoring. This device boasts a large LED-illuminated LCD screen, enabling users to view current, high, and low temperatures to the nearest tenth of a degree effortlessly. Measuring 3.5 inches in height, 7.25 inches in width, and 3.25 inches in depth, its compact design ensures flexibility in installation across various refrigeration units. The data logger is equipped with numerous features, such as temperature alarms with programmable tolerance bands to prevent unnecessary alerts, and an alarm delay function for user convenience. Its internal memory can store up to 32,000 data groups, and the included USB flash drive facilitates seamless data transfer to a PC or Mac. With reliable logging accuracy of $\pm 0.09^{\circ}\text{C}$ and a protected temperature sensor enclosed in glycol, this device ensures consistent performance. NIST calibrated and housed in durable steel, it promises both accuracy and longevity.

[Read More](#)

SKU: DL2BKIT



ECOLAB OR SOLUTIONS ORS-1075HS/ORS-2038D HUSH SLUSH SYSTEM

The **Ecolab ORS-1058HS hush slush machine** is designed to deliver velvet-soft slush essential for procedures requiring a delicate touch. Utilizing proprietary Auto-Stir technology, it automatically stirs the solution while cooling, eliminating the risk of sharp edges in conventional slush that could potentially cause trauma to delicate tissues. The machine features an Auto-Stir device that ensures consistent softness and a Disc-Drape that snaps onto the Auto-Stir device, providing optimal sterile coverage. Its caster wheels facilitate easy transportation and can be locked to minimize inadvertent movement, enhancing operational efficiency. The quiet operation of the ORS-1058HS ensures minimal disruption to the operating room staff. Included with the machine is the EcoLab OR Solutions Side by Side Warmer, ensuring comprehensive support during surgical procedures.

[Read More](#)

SKU: AM-ORS-1058HS



ACCUCOLD 24" WIDE ALL-FREEZER UP TO -30 DEG C CAPABLE

The **Accucold 24" Wide All-Freezer** is an advanced freezing solution designed to meet the stringent requirements of both commercial and medical facilities. With its capability to reach temperatures as low as -30°C (-22°F), this unit is particularly ideal for storing vaccines, pharmaceuticals, and preserving frozen foods over extended periods. It features a compact footprint of 23.63 inches in width and a height of 33.5 inches, making it suitable for installation in space-constrained areas. This freezer is built-in capable, thanks to its front-breathing system, and meets NSF-7 commercial standards, ensuring it adheres to the necessary safety and sanitation requirements. The freezer includes three slide-out drawers that help maintain a stable temperature environment each time the door is opened. Additional features such as a factory-installed keyed lock, a digital thermostat for precise temperature control, and an external temperature display with a programmable alarm enhance its functionality.



and security. The manual defrost system ensures consistent temperatures without fluctuations, promoting energy efficiency. Moreover, with its environmentally friendly design, using R134a refrigerant and being 100% CFC free, the Accucold 24" All-Freezer represents a reliable and sustainable choice for professional environments.

[Read More](#)

SKU: AM-VT65



ACCUCOLD 3.08 CU.FT. -85°C ULTRA LOW CHEST FREEZER 36.25IN H X 23.75IN W X 28.0IN D

Manufactured for Accucold in Denmark by ELCOLD, model LAB8 offers 3.08 cu.ft. of storage. Extra thick 4 1/2" insulation and static manual defrost operation help to maintain ultra low temperatures suitable for maintaining the integrity of biological samples, industrial materials, laboratory specimens, select vaccines, and other materials with highly specific cold storage requirements. Two inner lids offer added protection by shielding contents when the outer lid is opened. The digital microprocessor controls are mounted in the front for convenient access. This unit has a temperature range that can be set between -50 and -85°C, with a display viewable in Celsius or Fahrenheit. The thermostat's high/low temperature alarm warns if the interior goes out of the set range. There is also an audible 9V battery-powered alarm in the rear of the unit that will sound in the event of power failure.

[Read More](#)

SKU: AM-LAB8



ACCUCOLD 18 CU.FT. UPRIGHT VACCINE REFRIGERATOR, CERTIFIED TO NSF/ANSI 456 VACCINE STORAGE STANDARD

The **Accucold 18 Cu.Ft. Upright Vaccine Refrigerator** is a state-of-the-art solution designed specifically for the storage of vaccines, medications, and pharmacy applications, ensuring compliance with CDC/VFC guidelines and certification to NSF/ANSI 456 Vaccine Storage Standard. This refrigerator features a robust temperature management system, capable of maintaining an interior temperature with a minimal variation of $\pm 1^{\circ}\text{C}$, crucial for preserving the integrity of sensitive vaccines. It comes equipped with a 2-year calibration certificate, antimicrobial silver ion-coated handle to minimize germ spread, and adjustable temperature settings between $+2$ to $+8^{\circ}\text{C}$. Enhanced with a suite of alarms for temperature deviation, open door, power failure, and sensor disruption, it ensures the utmost safety and reliability. The unit also features remote alarm contacts for external monitoring, automatic defrost with adjustable cycles to prevent ice build-up, and a user-friendly microprocessor control panel with digital readouts for precise temperature management. The refrigerator offers ample storage space with adjustable plastic-coated wire shelves for optimal circulation, a factory-installed lock for security, and locking casters for mobility. Its efficient design, including a self-closing door, LED interior lighting, and a reversible door option, underscores its commitment to convenience and versatility in professional healthcare environments.

[Read More](#)

SKU: AM-ARS18PV456



ACCUCOLD UNDER COUNTER 6 CU.FT. VACCINE REFRIGERATOR, ADA HEIGHT

Accucold Under Counter 6 Cu.Ft. Vaccine Refrigerator, ADA Height

[Read More](#)

SKU: ARS62PVBIADA





24" WIDE ALL-REFRIGERATOR/ALL-FREEZER COMBINATION

Accucold offers innovative solutions to vaccine storage with stackable purpose-built refrigerator-freezer combinations design to support meeting CDC/VFC guidelines.

Sized just under 24" wide, the ARS12PV-FS24LSTACKMED2 combines the 12 cu.ft. ARS12PV auto defrost all-refrigerator with the compact FS24LMED2 all-freezer into one slim footprint using Accucold's StackRack. This ships in three boxes and requires some minor assembly upon delivery.

[Read More](#)

SKU: ARS12PV-FS24LSTACKMED2



ACCUCOLD 160 LB ICE & WATER DISPENSER WITHOUT FILTER

Experience the ultimate in convenience with the commercially approved Accucold countertop water and ice dispenser.

The AIWD160 is ETL-S listed to NSF-12 commercial standards. Measuring just 17.5 inches from the front leg to the unit rear makes this a perfect fit on standard counters, with a slim width under 14 inches. This unit can produce up to 160 pounds of ice per day. Three dispenser functions - ice, water, or ice+water - provide expanded functionality, while the illuminated display ensures easy access and control of the unit.

[Read More](#)

SKU: AIWD160



ACCUCOLD 160 LB ICE & WATER DISPENSER WITH FILTER

Experience the ultimate in convenience with the commercially approved Accucold countertop water and ice dispenser.

The AIWD160 is ETL-S listed to NSF-12 commercial standards. Measuring just 17.5 inches from the front leg to the unit rear makes this a perfect fit on standard counters, with a slim width under 14 inches. This unit can produce up to 160 pounds of ice per day. Three dispenser functions - ice, water, or ice+water - provide expanded functionality, while the illuminated display ensures easy access and control of the unit.

[Read More](#)

SKU: AIWD160FLTR



ACCUCOLD 450 LB ICE & WATER DISPENSER WITHOUT FILTER

Experience the ultimate in convenience with the commercially approved Accucold countertop water and ice dispenser.

The AIWD450 is ETL-S listed to NSF-12 commercial standards. Measuring just 20 inches from the front leg to the unit rear makes this a perfect fit on standard counters, with a slim 16.5 inch width. This unit can produce up to 378 pounds of ice per day. Three dispenser functions - ice, water, or ice+water - provide expanded functionality, while the illuminated display ensures easy access and control of the unit.

[Read More](#)

SKU: AIWD450



ACCUCOLD 450 LB ICE & WATER DISPENSER WITH FILTER

Experience the ultimate in convenience with the commercially approved Accucold countertop water and ice dispenser.

The AIWD450FLTR is ETL-S listed to NSF-12 commercial standards. Measuring just 20 inches from the front leg to the unit rear makes this a perfect fit on standard counters, with a slim 16.5 inch width. This unit can produce up to 378 pounds of ice per day. Three dispenser functions - ice, water, or ice+water - provide expanded functionality, while the illuminated display ensures easy access and control of the unit.

[Read More](#)

SKU: AIWD450FLTR



ACCUCOLD 4.8 CU.FT. -85°C ULTRA LOW CHEST FREEZER 36.25IN H X 28.5IN W X 28.0IN D

Accucold introduces a new series of ultra low temperature chest freezers designed for reliable storage at -85°C for laboratory, life science, and other institutional settings.

Manufactured for Accucold in Denmark by ELCOLD, model LAB11 offers 4.8 cu.ft. of storage. Extra thick 4 1/2" insulation and static manual defrost operation help to maintain ultra low temperatures suitable for maintaining the integrity of biological samples, industrial materials, laboratory specimens, select vaccines, and other materials with highly specific cold storage requirements. Two inner lids offer added protection by shielding contents when the outer lid is opened. The digital microprocessor controls are mounted in the front for convenient access. This unit has a temperature range that can be set between -50 and -85°C, with a display viewable in Celsius or Fahrenheit. The thermostat's high/low temperature alarm warns if the interior goes out of the set range. There is also an audible 9V battery-powered alarm in the rear of the unit that will sound in the event of power failure.

[Read More](#)

SKU: AM-LAB11



ACCUCOLD 4.8 CU.FT. -85°C ULTRA LOW CHEST FREEZER 36.25IN H X 28.5IN W X 28.0IN D WITH 9 STAINLESS VERTICAL BOX SHELVES

Accucold introduces a new series of ultra low temperature chest freezers designed for reliable storage at -85°C for laboratory, life science, and other institutional settings.

Manufactured for Accucold in Denmark by ELCOLD, model LAB11R9 offers 4.8 cu.ft. of storage. Extra thick 4 1/2" insulation and static manual defrost operation help to maintain ultra low temperatures suitable for maintaining the integrity of biological samples, industrial materials, laboratory specimens, select vaccines, and other materials with highly specific cold storage requirements. Two inner lids offer added protection by shielding contents when the outer lid is opened. The digital microprocessor controls are mounted in the front for convenient access. This unit has a temperature range that can be set between -50 and -85°C, with a display viewable in Celsius or Fahrenheit. The thermostat's high/low temperature alarm warns if the interior goes out of the set range. There is also an audible 9V battery-powered alarm in the rear of the unit that will sound in the event of power failure.

[Read More](#)

SKU: AM-LAB11R9



12 LB. DRAIN-FREE COUNTERTOP ICEMAKER CABINET STYLE

Sized at 15" wide, the BIM25 can produce up to 12 pounds of crescent-shaped ice per day, storing it in a removable insulated bin. This unit features a fully finished black cabinet with a stainless steel door for an attractive look. The front-breathing system allows the unit to be installed under counters. A side defrost drain is included for easier cleaning. No drain or pump is required for use. This unit includes a connection hose, but other hardware must be supplied by your plumber.

[Read More](#)

SKU: BIM25



SUMMIT 24" WIDE ALL-REFRIGERATOR, ADA COMPLIANT

Summit 24" Wide All-Refrigerator, ADA Compliant

[Read More](#)

SKU: AM-AL750W



ACCUCOLD ICE & WATER DISPENSER WITH FILTER

The AIWD282FLTR is ETL-S listed to NSF-12 commercial standards. Measuring 18 3/4 inches from the front leg to the unit rear makes this a perfect fit on standard counters. This unit can produce up to 282 pounds of ice per day. Three dispenser functions - ice, water, or ice+water - provide expanded functionality, while the illuminated display ensures easy access and control of the unit.

[Read More](#)

SKU: AIWD282FLTR




ACCUCOLD UNDER COUNTER MEDICATION/VACCINE REFRIGERATOR (GLASS DOOR)

Designed and purpose-built for pharmacy, medication, and vaccine applications to support meeting CDC/VFC guidelines, with antimicrobial handle

Exterior Dims: 32.5" H x 23.38" W x 24.38" D

[Read More](#)

SKU: AM-ARG6PV

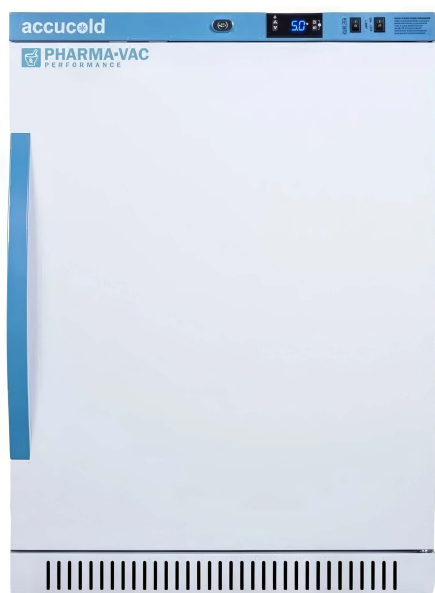

ACCUCOLD UNDER COUNTER MEDICATION/VACCINE REFRIGERATOR

Designed and purpose-built for pharmacy, medication, and vaccine applications to support meeting CDC/VFC guidelines, with antimicrobial handle

Exterior Dims: 32.5" H x 23.38" W x 24.38" D

[Read More](#)

SKU: AM-ARS6PV



ACCUCOLD 49 CU.FT. UPRIGHT PHARMACY REFRIGERATOR

The **Accucold 49 Cu.Ft. Upright Pharmacy Refrigerator** is a state-of-the-art refrigeration solution designed specifically for the stringent storage requirements of medical and pharmaceutical environments. With a substantial 49 cubic feet of storage capacity, this refrigerator offers ample space to accommodate a wide range of medical supplies, vaccines, and other temperature-sensitive materials. Engineered with precision, it features advanced temperature control systems to ensure a consistent and reliable environment, maintaining temperatures within the critical range necessary for preserving the efficacy of pharmaceuticals. The unit is equipped with high-quality shelving, digital temperature displays, and alarm systems to alert users of any deviations from the set parameters, thereby ensuring the utmost safety and compliance with regulatory standards. Its robust construction and energy-efficient design make the Accucold 49 Cu.Ft. Upright Pharmacy Refrigerator an indispensable asset for hospitals, clinics, and research facilities committed to maintaining the highest standards of medical storage.

[Read More](#)

SKU: AM-ARS49ML



ACCUCOLD 24" WIDE STACKED COMBINATION OFFERS SEPARATE REFRIGERATOR/FREEZER STORAGE

Accucold offers innovative solutions to vaccine storage with the STACKMED2 stackable refrigerator-freezer combination series.

Sized just under 24" wide, the FF7LW-VT65MLSTACKMED2 combines our popular FF7LWBIMED2 auto defrost all-refrigerator with the VT65MLBIMED2 all-freezer in one slim footprint. Made in Europe, the FF7LWBIMED2 offers 5.5 cu.ft. of storage, with fan-forced cooling and adjustable wire shelves to ensure an even temperature distribution. Model VT65MLBIMED2 is a low temperature all-freezer complete with interior drawers to ensure stored items are protected from ambient temperature when the door is open.

[Read More](#)

SKU: AM-FF7LW-VT65MLSTACKMED2



ACCUCOLD 49 CU.FT. UPRIGHT PHARMACY FREEZER

The **Accucold 49 Cu.Ft. Upright Pharmacy Freezer** is a state-of-the-art refrigeration solution designed to meet the stringent storage requirements of medical and pharmaceutical environments. With its substantial 49 cubic feet capacity, this freezer offers ample space for storing vaccines, medications, and other temperature-sensitive materials. Engineered with precision, it features advanced temperature control systems that ensure a consistent and reliable environment, maintaining temperatures as low as -25°C. The unit is equipped with high/low temperature alarms, a digital thermostat, and a lockable door to enhance security and monitoring. Its robust construction and energy-efficient design make it an ideal choice for hospitals, clinics, and research facilities that demand uncompromising performance and reliability in their cold storage solutions.

[Read More](#)

SKU: AM-AFS49ML

INDEX

1

- 1 Cu.Ft. Compact Medication/Vaccine Refrigerator, Certified to NSF/ANSI 456 Vaccine Storage Tabletop 3
- 12 lb. Drain-Free Countertop Ice maker Cabinet Style 11

2

- 24" Wide All-Refrigerator/All-Freezer Combination 8

A

- Accucold 160 lb Ice & Water Dispenser with Filter 9
- Accucold 160 lb Ice & Water Dispenser without Filter 8
- Accucold 18 Cu.Ft. Upright Vaccine Refrigerator, Certified to NSF/ANSI 456 Vaccine Storage Standard 6
- Accucold 220 lb Ice & Water Dispenser with Filter 2
- Accucold 24" Wide All-Freezer up to -30 deg C capable 4
- Accucold 24" Wide Stacked combination offers separate refrigerator/freezer storage 14
- Accucold 3.08 Cu.Ft. -85°C Ultra Low Chest Freezer 36.25in H x 23.75in W x 28.0in D 5
- Accucold 4.8 Cu.Ft. -85°C Ultra Low Chest Freezer 36.25in H x 28.5in W x 28.0in D 10
- Accucold 4.8 Cu.Ft. -85°C Ultra Low Chest Freezer 36.25in H x 28.5in W x 28.0in D with 9 Stainless Vertical Box Shelves 11
- Accucold 450 lb Ice & Water Dispenser with Filter 10
- Accucold 450 lb Ice & Water Dispenser without Filter 9
- Accucold 49 Cu.Ft. Upright Pharmacy Freezer 16
- Accucold 49 Cu.Ft. Upright Pharmacy Refrigerator 14
- Accucold Ice & Water Dispenser with Filter 12
- Accucold Product USB Digital Data Logger (DDL) 3
- Accucold Under Counter 6 Cu.Ft. Vaccine Refrigerator, ADA Height 7
- Accucold Under Counter Medication/Vaccine Refrigerator 13
- Accucold Under Counter Medication/Vaccine Refrigerator (Glass Door) 13

E

- ECOLAB OR SOLUTIONS ORS-1075HS/ORS-2038D HUSH SLUSH SYSTEM 4

S

- Summit 24" Wide All-Refrigerator, ADA Compliant 12