

TABLE OF CONTENTS

Scrub Sinks 2

 SS Series Dual Bay Surgical Scrub Sink with Infrared Operated Water & Soap 3

 Dual Bay Economy Sink with Manual Knee Operated Water 3

 Triple Bay Surgical Scrub Sink w/ Knee Operated Water & Soap 4

 Steris Manual Dual Bay Flexmatic Scrub Sink 4

 Dual Bay Flexematic Scrubsink, IR 4

 ES Series Scrub Sink Single Bay 5

 SS Series Dual Bay Surgical Scrub Sink W/ Knee Operated Water & Soap 5

 ES Series Dual Bay Economy Sink W/ Infrared Operated Water 6

 SS Series Single Bay Surgical Scrub Sink With Knee Operated Water & Soap 6

 Steris Amsco Flexmatic Single Bay Scrub Station 7



SCRUB SINKS

Built with safety in mind, these state-of-the-art surgical scrub sinks offer ease, convenience, and utility for all users. Durable and hands-free, Auxo Medical offers the highest quality scrub sinks from leading brands such as [Steris](#) and Mac Medical.

Offering one-, two-, and three-sink options, Auxo Medical has the right scrub sink for any operating room. Whether your facility needs a scrub sink with deep basins, adjustable spouts, easy to use mixing valves, or optimum water temperature and pressure, our surgical scrub sinks can fit your need. Our refurbished, stainless steel, easy to clean sinks offer user-friendly water controls with convenient knee and foot pedals. Auxo Medical provides only the most popular brands that focus on craftsmanship and industry-leading components.

Our catalog consists of brand new, reconditioned and refurbished equipment. Would you like to filter the equipment by:

[**BRAND NEW**](#)

[**RECONDITIONED**](#)

REFURBISHED



SS SERIES DUAL BAY SURGICAL SCRUB SINK WITH INFRARED OPERATED WATER & SOAP

The **SS Series Dual Bay Surgical Scrub Sink** is equipped with **infrared-operated water and soap dispensing systems**. Standard features include a deep sloping basin and an adjustable thermostatic controlled pressure regulating mixing valve for each station. Maintenance is simplified with check stop strainers. The sink features a high mount swivel gooseneck spout with an adjustable swivel spray aerator and a flow control valve. Backflow preventers are included to ensure water safety. Splash screen dividers are incorporated in double and triple station sink models. The sink is supported by a main support and Z-brackets, with optional in-wall support available. Constructed from 304 stainless steel, the sink is polished and ground to a #4 finish.

[Read More](#)

SKU: AM-SS64-IR




DUAL BAY ECONOMY SINK WITH MANUAL KNEE OPERATED WATER

The **Dual Bay Economy Sink with Manual Knee Operated Water**, part of the ES Series, offers a practical and hygienic handwashing solution. The ES47 model features a 47-inch wide sink with a deep, sloping basin designed for effective scrubbing. This sink is equipped with adjustable thermostatic controlled pressure regulating mixing valves (one per station) located beneath the sink, ensuring consistent water temperature. A high-mount swivel gooseneck spout with an adjustable swivel spray aerator provides comfortable water flow. It also includes a flow control valve and backflow preventers for efficient water management and safety. Constructed from durable 304 stainless steel with a #4 finish, the sink is easy to clean and maintain. Z-brackets are included for secure wall mounting (refer to specifications for in-wall support details). This model utilizes a manual knee-operated system for water activation, promoting hands-free operation and minimizing the risk of contamination.

[Read More](#)

SKU: AM-ES47

TRIPLE BAY SURGICAL SCRUB SINK W/ KNEE OPERATED WATER & SOAP

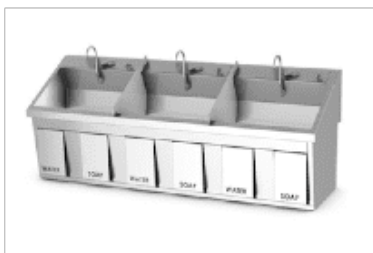

Triple Bay Surgical Scrub Sink w/ Knee Operated Water & Soap

[Read More](#)

SKU: AM-SS96



STERIS MANUAL DUAL BAY FLEXMATIC SCRUB SINK

Dual bay Flexmatic scrub station without timer or water control sensor. Each unit is constructed of welded 16 gauge stainless steel with a matte finish for easy cleaning and low maintenance. The built-in thermostatically controlled mixing valve provides optimum water temperatures up to 115° F (46° C)

[Download Brochure](#)

[Read More](#)

SKU: AM-FLEX-DB-M

DUAL BAY FLEXEMATIC SCRUBSINK, IR



Dual bay Flexmatic scrub station without timer or water control sensor. Each unit is constructed of welded 16 gauge stainless steel with a matte finish for easy cleaning and low maintenance. The built-in thermostatically controlled mixing valve provides optimum water temperatures up to 115° F (46° C). Chair carrier Included

Foot pedal or knee activated soap dispensers must be ordered separately (not included)

[Download Brochure](#)

[Read More](#)

SKU: AM-FLEX-DB-IR



ES SERIES SCRUB SINK SINGLE BAY

The ES Series Scrub Sink Single Bay is designed to address space and budget limitations without compromising functionality. Featuring a compact 25" wide design (ES25) with a deep, sloping basin, this scrub sink allows for effective and thorough handwashing. Standard features include an adjustable thermostatic controlled pressure regulating mixing valve mounted under the sink, a high-mount swivel gooseneck spout with an adjustable swivel spray aerator, a flow control valve, and backflow preventers. Infrared models are equipped with a 24-volt plug-in transformer. Constructed from durable 304 stainless steel polished and ground to a #4 finish, the ES Series scrub sink is both hygienic and aesthetically pleasing. Z-brackets are included for secure mounting (see specifications for in-wall support requirements).

[Download Brochure](#)

[Read More](#)

SKU: AM-ES25




SS SERIES DUAL BAY SURGICAL SCRUB SINK W/ KNEE OPERATED WATER & SOAP

The SS-Series Surgical Scrub Sinks were developed to answer space and economical constraints where needed. These sinks have smaller outside dimensions but still maintain a deep basin to allow for proper scrubbing.

[Download Brochure](#)

[Read More](#)

SKU: AM-SS64

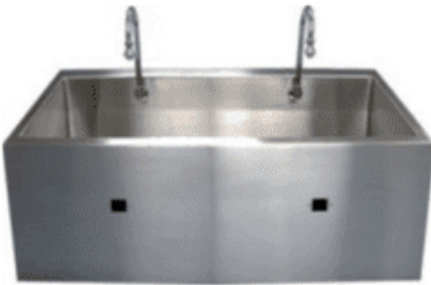

ES SERIES DUAL BAY ECONOMY SINK W/ INFRARED OPERATED WATER

The ES-Series Surgical Scrub Sinks were developed to answer space and economical constraints where needed. These sinks have smaller outside dimensions but still maintain a deep basin to allow for proper scrubbing.

[Download Brochure](#)

[Read More](#)

SKU: AM-ES25-IR



SS SERIES SINGLE BAY SURGICAL SCRUB SINK WITH KNEE OPERATED WATER & SOAP

The **SS Series Single Bay Surgical Scrub Sink** offers a blend of quality and affordability, making it an excellent choice for upgrading or equipping surgical facilities. This sink features knee-operated water and soap dispensing, promoting hygienic, hands-free operation. Standard features include a deep, sloping basin for effective scrubbing, an adjustable thermostatic controlled pressure regulating mixing valve per station for precise temperature control, and check stop strainers for easy maintenance. A high-mount swivel gooseneck spout with an adjustable swivel spray aerator ensures comfortable water flow. The sink also includes a flow control valve and backflow preventers for efficient water management. Constructed from 304 stainless steel with a polished and ground #4 finish, it provides durability and ease of cleaning. Splash screen dividers are included with double and triple station sink models. Main support and Z-brackets are provided for secure mounting, with optional in-wall support available.

[Download Brochure](#)

[Read More](#)

SKU: AM-SS32



STERIS AMSCO FLEXMATIC SINGLE BAY SCRUB STATION

Designed in accordance with Underwriters Laboratories standards to meet all of your scrub sink needs, Amsco® Flexmatic® Scrub Sinks have stood the test of time. These durable sinks feature continuous welded compartments and coved corners constructed entirely of stainless steel for easy cleaning. Durable, thermoformed Streamline scrub stations are made with built-in Microban antimicrobial protection with a seamless easy- to-clean basin. Easy-to-use knee panels make water control quick and convenient. The built-in thermostatically controlled mixing valve provides water temperatures up to 115° F (46° C). The standard "rose" spray head produces a shower-head-like stream of water without any aeration. A unique eyewash station attachment provides convenient operating room access to an emergency eye wash station. When water is turned on with knee panel, one pull of the center knob releases both orange safety caps to provide instant eye flushing. When knob is not activated, faucet functions as a scrub sink.

[Download Brochure](#)

[Read More](#)

SKU: AM-FLEX-SB-M

INDEX

D

Dual Bay Economy Sink with Manual Knee Operated Water 2

Dual Bay Flexematic Scrubsink, IR 4

E

ES Series Dual Bay Economy Sink W/ Infrared Operated Water 5

ES Series Scrub Sink Single Bay 4

S

SS Series Dual Bay Surgical Scrub Sink W/ Knee Operated Water & Soap 5

SS Series Dual Bay Surgical Scrub Sink with Infrared Operated Water & Soap 2

SS Series Single Bay Surgical Scrub Sink With Knee Operated Water & Soap 6

Steris Amsco Flexmatic Single Bay Scrub Station 6

Steris Manual Dual Bay Flexmatic Scrub Sink 3

T

Triple Bay Surgical Scrub Sink w/ Knee Operated Water & Soap 3