

TABLE OF CONTENTS

Sterile Processing Supplies & Accessories	3
SteriData Logger Paperless Sterilizer Cycle Printer For Statim	4
SteriData Logger Sterilizer Printer, Autoclave Printer, Paperless - Ritter / Midmark M9 / M11 Gen1 Series	4
Flush Tech Endoscope and lumened instrument flushing system	5
Tuttnauer TIVA Washer DIN Basket A002 Single	6
Unger Easy Adapter Fill Hose	6
SugiSonic Case of 12 replacement filters	6
Surgisonic 1220 (1) Verification test kit	7
Rennco LS24D Lift Seal Heat Sealer 110/115 volts	7
Printer Paper (1 Roll) for Steris Washer Disinfectors (P764330673)	8
Rennco LS24D Lift Seal Heat Sealer	8
Carbon Filter Refill Bags - Tuttnauer DS1000/DS1G Steam Distiller	9
SteriData Logger Sterilizer Printer, Autoclave Printer, Paperless - Ritter / Midmark M9 / M11 NEW Series	10
Pouch Sterilization Rack, Stainless Steel, Loop Style, 21 Slots 0.7in ea	10
Pouch Sterilization Rack, Stainless Steel, Loop Style, 8 Slots 2.5in ea	10
Pouch Sterilization Rack, Stainless Steel, Loop Style, 10 Slots 1.5in ea	11
Pouch Sterilization Rack, Stainless Steel, Loop Style, 5 Slots 1.5in ea	11
Qty 2: Black Cartridge Printer Ribbon For Steris Steam Sterilizers and Washers (P150828440)	12
Printer Paper for Steris Century & Evolution Sterilizers - 5 rolls/box (P387352558)	12
3M™ Attest™ Rapid Readout Biological Indicator 1291 for Steam Sterilization, 50/BX	13
Sharps Box, Combo package	14
Black Cartridge Printer Ribbon For Steris Steam Sterilizers and Washers (P150828440)	15
Printer Paper (1 Roll) for Steris Century and Evolution Sterilizers (P387352558)	15
Mucapur EnzymePower Cleaning Agent	15
Rennco LS18D Lift Seal Heat Sealer	16
Water-Mizer	17

Door Gasket for Castle 3433	17
Particulate and Softener 10in Filter System	17
Printer Ribbon Printon	18
Steris Dart Test (Daily Air Removal Test) - Single Unit	18
Thermal Printer Paper for Steris Sterilizers and Washers - Box of 5 (P129359008)	19
ESPE Attest Biological Indicators	19
ESPE Attest Biological Indicators by 3M-400/CS	20
3M Attest Auto Reader 390	20
Attest Steam Indicator Incubator	21
Healthmark TOSI Instrument Washer Test	22
Liquid Descaler	22
Instrument Lubricant 5GAL	22
Instrument Lubricant 4X1GAL	23
Auto Washer Detergent 4X1GAL	23
SERC Enzyme Detergent 5GAL	24
Nozzel, Injector, 444 Detergent	24
Dr. Weigert - Neodisher Z - Neutralizing Agent, Liquid Concentrate	25



STERILE PROCESSING SUPPLIES & ACCESSORIES

Auxo Medical offers a comprehensive range of sterile processing supplies and accessories designed to meet the diverse needs of medical and lab facilities. Our product lineup includes essential items such as heat sealers, printer ribbons, biological indicators, and specialized detergents, all crafted to ensure the highest standards of sterilization and cleanliness.

The Rennco LS24D Lift Seal Heat Sealer, for example, is a top-tier device that ensures a secure seal after a brief warm-up period, enhancing the efficiency and reliability of your sterilization process. Similarly, our 3M™ Attest™ Rapid Readout Biological Indicators provide a rapid and reliable method for monitoring steam sterilization, ensuring patient and personnel safety.

Auxo Medical not only focuses on providing high-quality products but also offers extensive support services to ensure the optimal performance of your sterile processing equipment. Our team is available for consultation, maintenance, and repair services across various locations, including Virginia, North Carolina, and several regions in Florida, among others. By choosing Auxo Medical, you are investing in a partner committed to enhancing the safety, efficiency, and reliability of your medical facility's sterile processing operations.

Our catalog consists of brand new, reconditioned and refurbished equipment. Would you like to filter the equipment by:

[BRAND NEW](#)

[RECONDITIONED](#)

[REFURBISHED](#)



STERIDATA LOGGER PAPERLESS STERILIZER CYCLE PRINTER FOR STATIM

SteriData™ Logger for SciCan Statim 2000 and 5000 cassette sterilizers*
(must have must have serial port in back and current software)

The paperless alternative to a paper Printer

This SteriData Logger unit plugs in just like the paper printer and sits on top of your Sterilizer. Now all your sterilizer cycle data is saved to the SD card (included) automatically.

When the cycle data is required, simply take the SD card to your computer where all the cycle files are arranged by year and month, date and time.

The individual sterilizer cycle files are recorded as text files and easily displayed and/or printed out from your computer as required. Optionally, you can include your practice name and the sterilizer serial number which will be added to every cycle file.

[Read More](#)

SKU: AM-S1219DL

Price: \$1,075.00



STERIDATA LOGGER STERILIZER PRINTER, AUTOCLAVE PRINTER, PAPERLESS - RITTER / MIDMARK M9 / M11 GEN1 SERIES

SteriData™ Logger for Midmark Ritter M9 and M11 Ultraclaves, M3 Ultrafast sterilizer, and SciCan Statim 2000 and 5000 cassette sterilizers*

The paperless alternative to a paper Printer

This SteriData Logger unit plugs in just like the paper printer and sits on top of your Sterilizer. Now all your sterilizer cycle data is saved to the SD card (included) automatically.

When the cycle data is required, simply take the SD card to your computer where all the cycle files are arranged by year and month, date and time.

The individual sterilizer cycle files are recorded as text files and easily displayed and/or printed out from your computer as required. Optionally, you can include your practice name and the sterilizer serial number which will be added to every cycle file.

[Read More](#)

SKU: AM-M1219DL-G1

Price: \$1,095.00



FLUSH TECH ENDOSCOPE AND LUMENED INSTRUMENT FLUSHING SYSTEM

The **FlushTech™ GI Endoscope and Lumened Instrument Flushing System (7260030)** is a compact, durable, dual-port device designed to improve the efficiency and compliance of endoscope reprocessing. Its robust system facilitates the thorough flushing of detergent, rinse water, and air through all internal and elevator wire channels of flexible endoscopes, reducing reprocessing time while ensuring adherence to strict flushing requirements. The unit features a user-friendly interface, a flow rate of 0.75 GPM, and a maximum suction height of 6 feet, with compatibility for industry-standard high-level disinfectants and detergents in diluted form—excluding isopropyl alcohol. Compact and lightweight at 9.4 lbs, with an IPX0 case rating, the FlushTech system maintains detailed processing reports to support quality assurance and survey audits, making it an essential tool for efficient and compliant endoscope reprocessing.

[Read More](#)

SKU: AM-FTGI



Tuttnauer TIVA Washer DIN Basket A002 Single

TUTTNAUER TIVA WASHER DIN BASKET A002 SINGLE

The **Tuttnauer TIVA Washer DIN Basket A002** is a durable, stainless steel cleaning basket designed for use in sterilization and washer systems. Measuring **19 inches x 10 inches x 2 inches**, this single basket provides a secure and organized container for small surgical instruments, trays, or other medical devices during cleaning cycles. Its DIN (Deutsches Institut für Normung) standards ensure compatibility and reliable performance within compatible washer/disinfector systems, helping facilitate thorough and efficient cleaning processes.

[Read More](#)

SKU: ATS238-0293



UNGER EASY ADAPTER FILL HOSE

The **Unger Easy Adapter Fill Hose** is a 6-foot long, durable hose designed to streamline water filling in cleaning operations. Its length provides easy reach from faucets to buckets, minimizing movement and effort for users, while its robust construction ensures long-lasting performance in high-traffic environments. The hose features a quick and secure faucet connection that prevents detachment, supporting continuous water flow for uninterrupted cleaning tasks. Compatible with Unger dual-bucket systems, it enhances water management efficiency, reduces the need for frequent faucet adjustments, and boosts overall productivity and convenience during cleaning procedures.

[Read More](#)

SKU: AM-WH-180



SugiSonic Case of 12 replacement filters

SUGISONIC CASE OF 12 REPLACEMENT FILTERS

The **SurgiSonic Case of 12 Replacement Filters** includes in-line 70-micron filters designed to be used with ultrasonic cleaning systems. These filters serve as essential filtering assemblies that immediately trap dislodged debris during cleaning, preventing contamination of the ultrasonic tank's water solution. Replacing these filters regularly ensures optimal cleaning performance, maintains system efficacy, and upholds sterilization standards by minimizing particulate buildup in the ultrasonic process.

[Read More](#)

SKU: AM-SS-3112



SURGISONIC 1220 (1) VERIFICATION TEST KIT

The **SurgiSonic 1220 Verification Test Kit** is an essential set designed to ensure the proper functioning and performance of ultrasonic cleaning equipment. It includes a vial bracket to secure test vials within the cleaning basket, five ultrasonic test vials for evaluating cleaning effectiveness, a temperature strip to monitor process temperature, a lumen test strip holder, and five lumen test strips to assess cleaning in lumened or cannulated instruments. Additionally, the kit comes with a record-keeping sheet for documentation and instructions for use (IFU). This comprehensive kit helps maintain compliance with sterilization standards and confirms optimal operation of ultrasonic cleaners.

[Read More](#)

SKU: AM-SS-1220



RENNCO LS24D LIFT SEAL HEAT SEALER 110/115 VOLTS

The **Rennco LS24D Lift Seal Heat Sealer (110/115V)** is a compact, easy-to-use 24" wide unit that delivers a consistent 1/2" seal with about 30 psi of sealing pressure in roughly 2 seconds. It features a digital temperature control, a safety lift bar with automatic lower bar release, and a non-stick heated sealing surface for reliable operation. Designed for 110/115 V, it draws about 5 A and includes a digital readout for precise settings. With a sturdy 60 lb body and simple edge-in sealing, this sealer is ideal for quick, repeatable bag sealing in small production or packaging environments.

[Read More](#)

SKU: AM-LS24D-115

Price: \$3,995.00



PRINTER PAPER (1 ROLL) FOR STERIS WASHER DISINFECTORS (P764330673)

This **printer paper roll** ensures that your equipment operates efficiently and maintains accurate records of the disinfection processes, contributing to the overall reliability and performance of your sterilization procedures.

Compatible for use in various STERIS Corporation washer disinfectors. Each roll measures 70 feet in length and is compatible with a range of models including the 1300 Cart Washer, 3052 Single Chamber Washer, 5052 Single Chamber Washer, 7053L Single Chamber Washer, Finn-Aqua, Reliance 1227, Reliance 300XLS, and the Reliance Vision Single Chamber Washer Disinfectors. **Steris OEM Part Number: 764330-673**

QTY: 1 Roll

[Read More](#)

SKU: 764330-673

Price: \$11.95



RENNCO LS24D LIFT SEAL HEAT SEALER

Rennco offers the 24 LSD heat sealer. This lift seal heat sealer is ready for use after a short "warm-up" period. You simply adjust the temperature control knob on the heat sealing equipment for the material being used, place the edge of the bag in the seal area, lift the safety bar, and the lift seal automatically makes a 1/2" seal with approximately 30 PSI of consistent sealing pressure. In approximately 2 seconds, the lower safety bar on the heat sealer automatically opens for easy removal of the bag. Constant temperature is maintained by an internal digital heat control during the sealing process. The lift seal heat sealer will not cycle if an object exceeding 1/8" thick is in the sealing area, providing safe operation.

The Rennco 24 LS lift seal heat sealer is quick as well as easy to operate and maintain. The lower silicone rubber bar is raised by a rugged mechanical drive to meet the stationary 1/2" wide heated sealing bar above. This bar is covered with a special "non-stick" fabric which is manually advanced by a simple turn of a knob as the fabric becomes worn. The fabric roll is ample for months of use and is quickly replaced by snapping the new roll into position.

[Read More](#)

SKU: AM-LS24D-120

Price: \$3,995.00



CARBON FILTER REFILL BAGS - TUTTNAUER DS1000/DS1G STEAM DISTILLER

These specialized refill bags contain activated carbon, which is a powerful adsorbent known for its ability to remove impurities, including chlorine, sediments, and volatile organic compounds (VOCs) from water. By integrating these carbon filters into the distillation process, users of the Tuttnauer DS1000 / DS1G systems can achieve higher quality distilled water, which is essential for laboratory applications and medical purposes. The convenience of the refill bags also means that maintaining the distiller's performance and water quality over time is straightforward, ensuring that users consistently have access to clean, pure distilled water without the hassle of frequent equipment maintenance.

Qty 6 - Compatible with DS1000 or DS1G

[Read More](#)

SKU: 8606



STERIDATA LOGGER STERILIZER PRINTER, AUTOCLAVE PRINTER, PAPERLESS - RITTER / MIDMARK M9 / M11 NEW SERIES

SteriData™ Logger for Midmark Ritter M9 and M11 Ultraclaves, M3 Ultrafast sterilizer, and SciCan Statim 2000 and 5000 cassette sterilizers*

The paperless alternative to a paper Printer

This SteriData Logger unit plugs in just like the paper printer and sits on top of your Sterilizer. Now all your sterilizer cycle data is saved to the SD card (included) automatically.

[Read More](#)

SKU: AM-M1219-DL

Price: \$1,095.00



**Pouch Sterilization Rack,
Stainless Steel, Loop Style,
21 Slots 0.7in ea**

POUCH STERILIZATION RACK, STAINLESS STEEL, LOOP STYLE, 21 SLOTS 0.7IN EA

21 slots
0.7in spacing
16.5in x 4in x 4in

[Read More](#)

SKU: AM-PR-L1744-21



POUCH STERILIZATION RACK, STAINLESS STEEL, LOOP STYLE, 8 SLOTS 2.5IN EA

**Pouch Sterilization Rack,
Stainless Steel, Loop Style,
8 Slots 2.5in ea**

8 slots
2.5in spacing
20.5in x 6in x 6in

[Read More](#)

SKU: AM-PR-L2066-8



POUCH STERILIZATION RACK, STAINLESS STEEL, LOOP STYLE, 10 SLOTS 1.5IN EA

**Pouch Sterilization Rack,
Stainless Steel, Loop Style,
10 Slots 1.5in ea**

10 slots
1.5in spacing
16.5in x 4in x 4in

[Read More](#)

SKU: AM-PR-L1844-10



POUCH STERILIZATION RACK, STAINLESS STEEL, LOOP STYLE, 5 SLOTS 1.5IN EA

**Pouch Sterilization Rack,
Stainless Steel, Loop Style,
5 Slots 1.5in ea**

5 slots
1.5in spacing
8.7in x 4in x 4in

[Read More](#)

SKU: AM-PR-L0944-5



QTY 2: BLACK CARTRIDGE PRINTER RIBBON FOR STERIS STEAM STERILIZERS AND WASHERS (P150828440)

Epson ERC-22 B Black Cartridge Printer ribbon For:

Steris Century Steam Sterilizers
Steris 400 Series Sterilizers
Steris Evolution Sterilizers
Steris Vision Washers
Steris 3052, 5052, 2532 washers
Steris LAB110, 250 Sterilizers
Steris 250LS 630 LS sterilizers

Steris Part Number: 150828-440 and P150828-440

[Read More](#)

SKU: 150828-440

Price: \$18.59



PRINTER PAPER FOR STERIS CENTURY & EVOLUTION STERILIZERS - 5 ROLLS/BOX (P387352558)

Printer Paper for use in Steris Century and Evolution Sterilizers. Sold in Box of 5 Rolls.

Fits Steris Century sterilizer Healthcare and Scientific Models: V116, SG116, V120, SG120, V136, SG136, V148, SG148, V160, SG160, SV1262
Fits Steris Evolution sterilizer Healthcare and Scientific Models: HC600, HC900, HC1200

Steris OEM Part#: P387352558

[Read More](#)

SKU: 387352-558

Price: \$24.95



3M™ ATTEST™ RAPID READOUT BIOLOGICAL INDICATOR 1291 FOR STEAM STERILIZATION, 50/BX

3M Attest 1291 Rapid Readout Biological Indicators For Gravity Steam Sterilizers

[Read More](#)

SKU: AM-MMM1291

Price: \$288.80



SHARPS BOX, COMBO PACKAGE

[Read More](#)

SKU: AM-PLA145002

Price: \$116.70





BLACK CARTRIDGE PRINTER RIBBON FOR STERIS STEAM STERILIZERS AND WASHERS (P150828440)

Epson ERC-22 B Black Cartridge Printer ribbon For:

Steris Century Steam Sterilizers
Steris 400 Series Sterilizers
Steris Evolution Sterilizers
Steris Vision Washers
Steris 3052, 5052, 2532 washers
Steris LAB110, 250 Sterilizers
Steris 250LS 630 LS sterilizers

Related Steris Part Number: 150828-440 and P150828-440 *

*Please note this listing is for a single (1) cartridge

[Read More](#)

SKU: AM-ERC22

Price: \$9.95



PRINTER PAPER (1 ROLL) FOR STERIS CENTURY AND EVOLUTION STERILIZERS (P387352558)

This **printer paper**, that is **compatible with Steris Century and Evolution Sterilizers**, ensures clear and accurate printouts of cycle data, which are essential for record-keeping and compliance with industry standards. By choosing our printer paper, you're investing in a dependable product that supports your sterilization documentation needs, helping you maintain a high standard of care and operational excellence.

QTY: 1 Roll

[Read More](#)

SKU: AM-704449

Price: \$4.99



MUCAPUR ENZYMEPOWER CLEANING AGENT

The Mucapur EnzymePower Cleaning Agent is fit for automated reprocessing, manual pre-cleaning, or soaking of sensitive surgical instruments. The product provides outstanding cleaning efficiency with the gentle power of bioactive enzymes.

Features:

- PH-neutral enzymatic cleaning agent
- High material compatibility
- 5-liter container

[Read More](#)

SKU: AM-94069027

Price: \$229.23



RENNCO LS18D LIFT SEAL HEAT SEALER


The Rennco Model 18 LS heat sealer is ready for use after a short "warm-up" period. Simply enter the temperature setting for the material being used, place the edge of the bag in the seal area, lift the safety bar, and the lift seal automatically makes a 1/2" seal with approximately 30 PSI of consistent sealing pressure. In approximately 2 seconds, the lower safety bar on the heat sealer automatically opens for easy removal of the bag. Constant temperature is maintained by an internal digital heat control during the sealing process. The lift seal heat sealer will not cycle if an object exceeding 1/8" thick is in the sealing area, providing safe operation.

The Rennco 18 LS Lift Seal heat sealer is quick as well as easy to operate and maintain. The lower silicone rubber bar is raised by a rugged mechanical drive to meet the stationary 1/2" wide heated sealing bar above. This bar is covered with a special "non-stick" fabric which is manually advanced by a simple turn of a knob as the fabric becomes worn. The fabric roll is ample for months of use and is quickly replaced by snapping the new roll into position.

[Read More](#)

SKU: AM-LS18D-120

Price: \$2,999.95

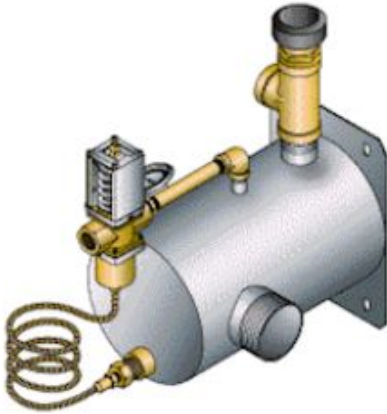

WATER-MIZER

The Water-Mizer, a water-saving device, eliminates the "tempering" water consumption during the non-sterilizing portion of a sterilizer's cycle, thereby substantially reducing overall water consumption.


Features:

- Shipping Dimensions: 20"Lx18"Wx15"H
- Weight: 25 lbs

[Read More](#)

SKU: AM-RPA466

Price: \$2,915.54

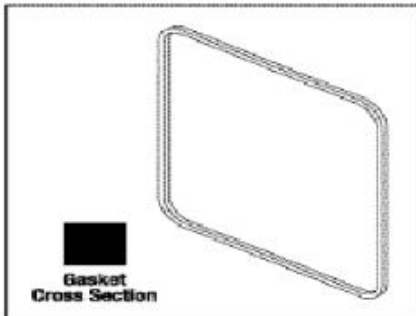

DOOR GASKET FOR CASTLE 3433

The Door Gasket for Castle 3433 significantly cuts down on the infiltration of sound, light, drafts, and moisture.

Features:

Shipping Dimensions: 30"Lx30"Wx1"H

Weight: 10lbs

[Read More](#)

SKU: 60864

Price: \$165.90



PARTICULATE AND SOFTENER 10IN FILTER SYSTEM

[Read More](#)

SKU: AM-WC-PS-10

Price: \$428.00

PRINTER RIBBON PRINTON

The Printer Ribbon Printon is fit for use in sterilizer printers.

Features:

- Printer Ribbon, Printon, Olivetti
- Shipping Dimensions: 6"Lx6"Wx2"H
- Weight: 1 lb

[Read More](#)

SKU: 129356-118

Price: \$54.92





STERIS DART TEST (DAILY AIR REMOVAL TEST) - SINGLE UNIT

The Dart Test Single Unit is compliant with ANSI/AAMI/ISO 11140-4 requirements for class 2 air removal indicator systems.

Features:

- Can operate at 270-274°F/132-134°C with an exposure time of 3.5 to 4 minutes
- No more spending time looking for towels and preparing packs
- Can eliminate all the waste associated with paper box products
- Light, compact, one-piece system arrives ready for use

[Read More](#)

SKU: NB215

Price: \$22.20



THERMAL PRINTER PAPER FOR STERIS STERILIZERS AND WASHERS - BOX OF 5 (P129359008)

Box of 5 Thermal Printer Paper for use in Steris Stage 3 sterilizer and Steris 444 and Synergy washer disinfectors.

Steris OEM Part#: P129359008

Sterilizer Models: 3011, 3013, 3021, 3023, 3031, 3033, 3041, 3043, 3051, 3053

Washer Models: 444, Synergy, Genfore, 430

- Shipping Dimensions: 10"Lx6"Wx2"H

- Weight: 2 lb

[Read More](#)

SKU: 129359-008

Price: \$19.95

ESPE ATTEST BIOLOGICAL INDICATORS



The Espe Attest Biological Indicators provides a maximum challenge to various steam sterilization cycles.

Features:

- This Biological Indicator is for use with steam at 270(deg)F/132(deg)C
- Gravity 270 only 100/bx
- Self-contained biological indicators (25/box)

[Read More](#)

SKU: AM-MMM1261Z

Price: \$301.93

ESPE ATTEST BIOLOGICAL INDICATORS BY 3M-400/CS



The Espe Attest Biological Indicators provides a maximum challenge to various steam sterilization cycles.

Features:

- This Biological Indicator is for use with steam at 270(deg)F/132(deg)C
- For Gravity 270 by 3M - 400/CS
- Self-contained biological indicators (25/box)
- Design provides a maximum challenge to various steam sterilization

[Read More](#)

SKU: AM-MMM1261

Price: \$1,201.56

3M ATTEST AUTO READER 390



The 3M Attest 390 Auto-reader for Steam reads the fluorescence produced by positive 3M Attest Rapid Readout Biological Indicators 1291 and 1292 to indicate a sterilization process failure. Featuring advanced electronics and a backlit LCD panel, the 3M Attest Auto-reader 390 turns up the dial on reliability and accuracy. It gives you even more confidence in your efforts to ensure patient safety while working with your existing 3M Attest Rapid Readout products.

Features:

- Includes ethernet cable for web app connectivity.
- Audible alarm sounds for improved monitoring
- Designed for use with Attest Rapid Readout Biological Indicators and Test Packs 1291, 1292, 1296, 1296F, 41382, 41382F
- Intelligent display with time-to-result countdown feature
- Display can be viewed remotely using your Smartphone or PC
- Manage record-keeping with web app and your internet connected PC
- Connects to select sterile processing software systems; ask your software system representative for details

[Read More](#)

SKU: AM-3M-AAR-390

Price: \$1,823.96

ATTEST STEAM INDICATOR INCUBATOR



The Attest Steam Indicator Incubator is a convenient, reliable system for monitoring steam sterilization.

Features:

- Test results are known in 24-48 hours
- Operates at 56(deg)C - 120 V.A.C.
- 14-vial capacity ideal for large and small offices.

[Read More](#)

SKU: AM-MMM116

Price: \$289.64

HEALTHMARK TOSI INSTRUMENT WASHER TEST



The Healthmark TOSI Instrument Washer Test is the first device to provide a consistent, repeatable, and reliable method for evaluating the cleaning effectiveness of the automated instrument washer.

Features:

- Easy-to-use, blood-soil device that helps correlate the cleaning challenge of surgical instruments
- Clear, plastic cover for viewing of test results
- Visually clean plate helps indicate effective cleaning

Sold in case of 30

[Read More](#)

SKU: AM-WT-101

Price: \$609.67

LIQUID DESCALER



The Liquid Descaler acts as an excellent chelating agent on iron sulfates and iron chlorides. dissolve carbonate, sulfate, magnesium, and iron deposits.

Features:

- Parts source # 105108 4/1 gallon SERC # 105108 4/1gallon
- Precision Dynamics # DESCALE 518 5gallon container
- Shipping Dimensions: 32"Lx20"Wx14"H
- Weight: 30 lbs

[Read More](#)

SKU: 105108

Price: \$269.40



INSTRUMENT LUBRICANT 5GAL

The Instrument Lubricant 5GAL is a non-silicone based lubricant with a rust inhibitor that prolongs instrument life. It is not for use as a medical device high-level disinfectant or sterilant. The product is phosphate-free.

Features:

- Shipping Dimensions: 24"Lx24"Wx24"H
- Weight: 40 lbs

[Read More](#)

SKU: S-301-5

Price: \$161.00



INSTRUMENT LUBRICANT 4X1GAL

The Instrument Lubricant 4X1GAL is a non-silicone based lubricant with a rust inhibitor that prolongs instrument life. It is not for use as a medical device high-level disinfectant or sterilant. The product is phosphate-free.

Features:

- Shipping Dimensions: 24"Lx24"Wx24"H
- Weight: 40 lbs

[Read More](#)

SKU: S-301-CG

Price: \$150.85

AUTO WASHER DETERGENT 4X1GAL



The Auto Washer Detergent 4X1GAL

[Read More](#)

SKU: S-400-CG

Price: \$112.89

SERC ENZYME DETERGENT 5GAL



The SERC Enzyme Detergent 5GAL catalyzes the reaction between stains and the water solution, thus aiding stain removal. View more products!

[Read More](#)

SKU: S-504-5

Price: \$257.87

NOZZEL, INJECTOR, 444 DETERGENT



The Nozzel, Injector, 444 detergent is a special part of the washing and disinfection process.

Shipping Dimensions: 6"Lx6"Wx4"H

Weight: 5 lbs

[Read More](#)

SKU: 117954-320

Price: \$38.28



DR. WEIGERT - NEODISHER Z - NEUTRALIZING AGENT, LIQUID CONCENTRATE

Dr. Weigert Neodisher Z is a specialized cleaning agent formulated for the effective removal of stubborn residues from laboratory glassware and other instruments. Known for its superior performance in eliminating inorganic and organic contaminants, Neodisher Z is widely used in scientific and medical settings to ensure the highest standards of cleanliness and sterility. Its unique composition is designed to target and dissolve residues without causing damage to delicate surfaces, making it an essential component in the maintenance of precise and reliable laboratory environments. Additionally, Neodisher Z is highly regarded for its ease of use and compatibility with various cleaning protocols, contributing to the overall efficiency and safety of laboratory operations.

5 liter container

[Read More](#)

SKU: AM-US483033

Price: \$229.95

INDEX

3

3M Attest Auto Reader 390 19

3M™ Attest™ Rapid Readout Biological Indicator 1291 for Steam Sterilization, 50/BX 11

A

Attest Steam Indicator Incubator 19

Auto Washer Detergent 4X1GAL 22

B

Black Cartridge Printer Ribbon For Steris Steam Sterilizers and Washers (P150828440) 13

C

Carbon Filter Refill Bags - Tuttnauer DS1000/DS1G Steam Distiller 7

D

Door Gasket for Castle 3433 15

Dr. Weigert - Neodisher Z - Neutralizing Agent, Liquid Concentrate 23

E

ESPE Attest Biological Indicators 18

ESPE Attest Biological Indicators by 3M-400/CS 18

F

Flush Tech Endoscope and lumened instrument flushing system 3

H

Healthmark TOSI Instrument Washer Test 20

I

Instrument Lubricant 4X1GAL 21

Instrument Lubricant 5GAL 21

L

Liquid Descaler 20

M

Mucapur EnzymePower Cleaning Agent 14

N

Nozzel, Injector, 444 Detergent 23

P

Particulate and Softener 10in Filter System 16

Pouch Sterilization Rack, Stainless Steel, Loop Style, 10 Slots 1.5in ea 9

Pouch Sterilization Rack, Stainless Steel, Loop Style, 21 Slots 0.7in ea 8

Pouch Sterilization Rack, Stainless Steel, Loop Style, 5 Slots 1.5in ea 10

Pouch Sterilization Rack, Stainless Steel, Loop Style, 8 Slots 2.5in ea 9

Printer Paper (1 Roll) for Steris Century and Evolution Sterilizers (P387352558) 13

Printer Paper (1 Roll) for Steris Washer Disinfectors (P764330673) 6

Printer Paper for Steris Century & Evolution Sterilizers - 5 rolls/box (P387352558) 11

Printer Ribbon Printon 16

Q

Qty 2: Black Cartridge Printer Ribbon For Steris Steam Sterilizers and Washers (P150828440) 10

R

Rennco LS18D Lift Seal Heat Sealer 14

Rennco LS24D Lift Seal Heat Sealer 7

Rennco LS24D Lift Seal Heat Sealer 110/115 volts 6

S

SERC Enzyme Detergent 5GAL 22

Sharps Box, Combo package 12

SteriData Logger Paperless Sterilizer Cycle Printer For Statim 2

SteriData Logger Sterilizer Printer, Autoclave Printer, Paperless - Ritter / Midmark M9 / M11 Gen1
Series 3

SteriData Logger Sterilizer Printer, Autoclave Printer, Paperless - Ritter / Midmark M9 / M11 NEW Series
8

Steris Dart Test (Daily Air Removal Test) - Single Unit 17

SugiSonic Case of 12 replacement filters 5

Surgisonic 1220 (1) Verification test kit 5

T

Thermal Printer Paper for Steris Sterilizers and Washers - Box of 5 (P129359008) 17

Tuttnauer TIVA Washer DIN Basket A002 Single 4

U

Unger Easy Adapter Fill Hose 4

W

Water-Mizer 15