

Q'STRAIT[®]

AMBULOK[®]

UNIVERSAL COT SECUREMENT SYSTEM


FOR **NON-EMERGENCY TRANSPORT**
IN **CARGO VEHICLES**

INSTANTLY BRING YOUR FLEET
INTO **STANDARDS COMPLIANCE**

Maximize Safety

Improve Productivity

Reduce Operating Costs

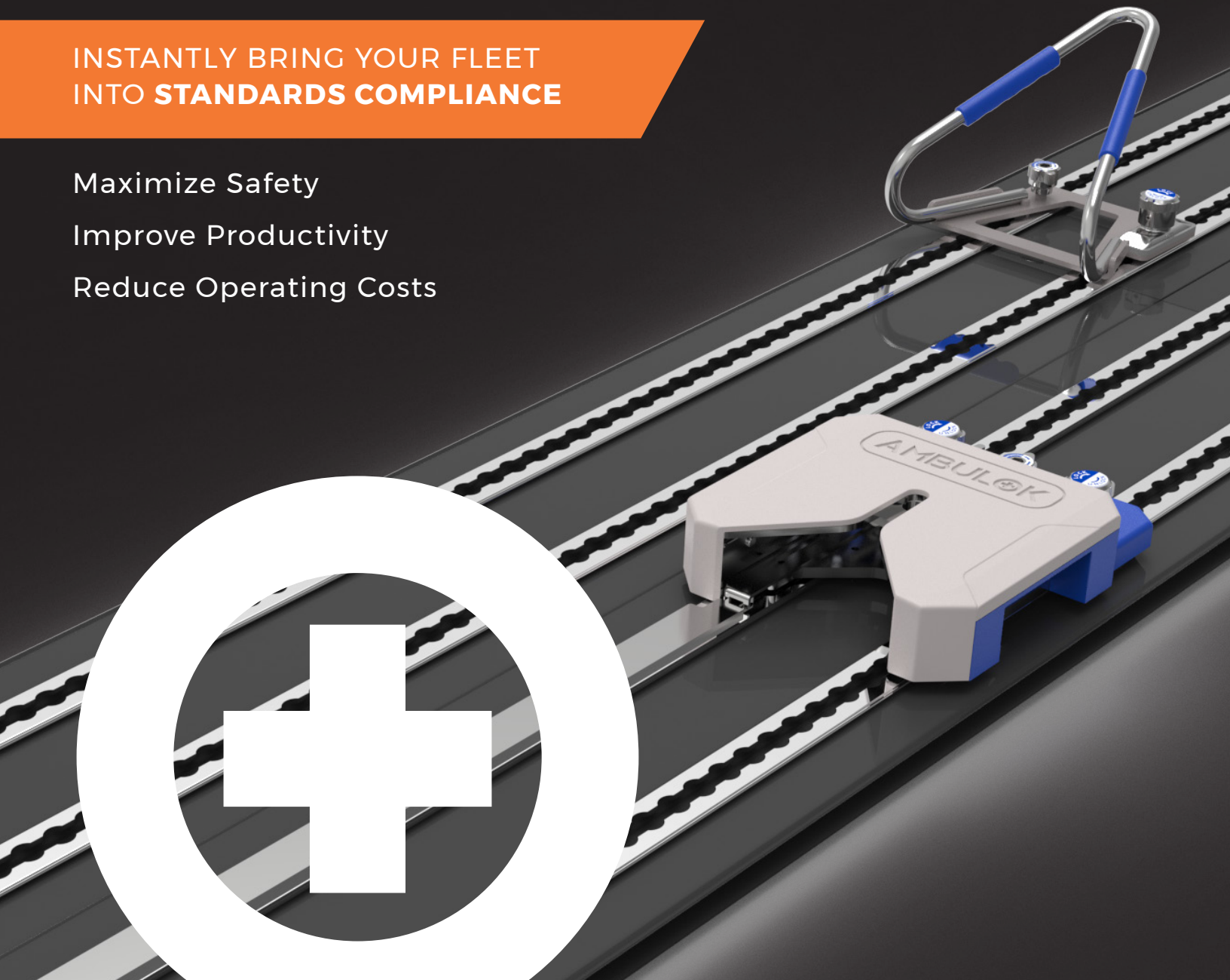

INTRODUCING

AMBULOK⁺ THE FUTURE OF NEMT COT SECUREMENT



INSTANT COMPLIANCE

The AMBULOK system can be installed into any existing NEMT vehicle with OMNI Floor. The innovative cot bracket system lets you **secure virtually any cot**, including those already in your fleet, instantly bringing that cot into standard compliance in the most efficient and economical way possible.

FULLY CRASH CERTIFIED

AMBULOK meets **SAE J3027**, the crash standard for ambulance-based equipment. Using proven technology adapted from our world-class wheelchair securement applications, AMBULOK safely secures cots in transport, protecting patients and operators from major injuries - and your business from liability claims.



EFFORTLESS FLEXIBILITY

The low profile AMBULOK docking system requires **no tools** for install or removal, allowing you to quickly transition between cot and wheelchair transport configurations as needed. The entire process **takes just minutes** and makes the cleaning, storage and upkeep of the system a breeze.



HIGHLY INTUITIVE

Operation of the AMBULOK system was patterned around modern First Responders, with each step designed to use 'medic muscle memory', or movements that both EMS and NEMT techs are already familiar with. Every part is lightweight (just 10 lbs. or less) and ergonomic, enabling easy **one-handed operation**.




Not only does the AMBULOK system bring **ANY Cot** into standard compliance, it also works seamlessly with **the most popular cots in the industry**.

THE COMPLETE AMBULOK SYSTEM INCLUDES:



DOCKING SYSTEM

Modular assembly that can be rapidly
installed and removed

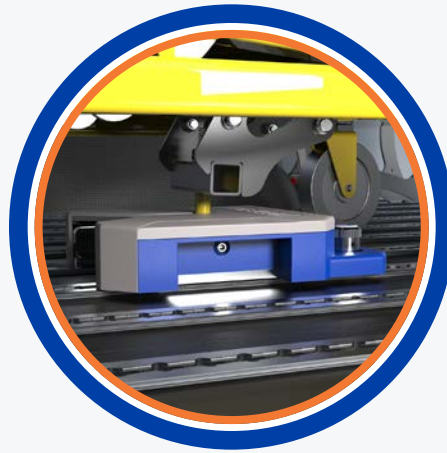






COT SECUREMENT BRACKET

Docking system interface that installs with
no cot modification needed



PATIENT RESTRAINT

4-point patient harness with one-hand
release buckle

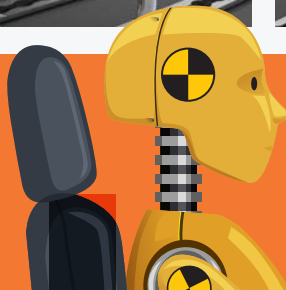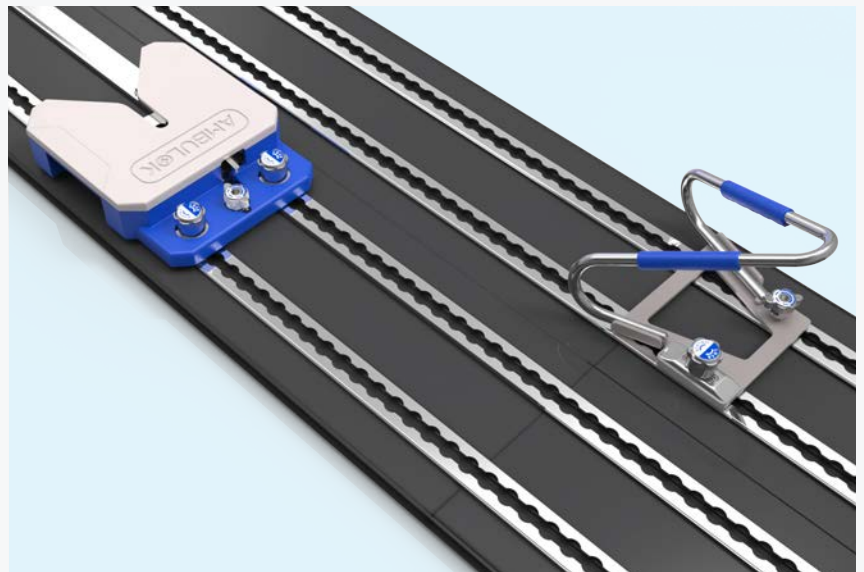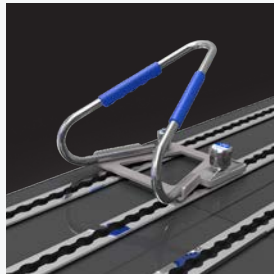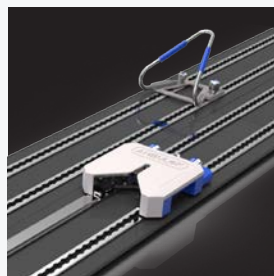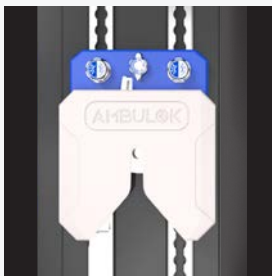

WHAT IS SAE J3027?

SAE J3027 recommends the testing procedures required to evaluate the integrity of a ground ambulance-based patient cot, cot securement system, and patient restraint when exposed to a frontal, side, or rear impact.

ASK ABOUT OTHER APPLICATIONS
**AMBULOK FOR
NON-EMERGENCY
MINIVAN VEHICLES**

Available for Rear-Entry Chrysler Voyager





AMBULOK CARGO SPECIFICATIONS

DOCKING SYSTEM SPECS

Length: 9.12" // Width: 9.12" // Height: 2.62" // Weight: 10 lbs.

RELEASE MECHANISM SPECS

Length: 35" // Width: 1.22" // Height: 2.62" // Weight: 1.25 lbs.

TOTAL COT SECUREMENT ZONE SPECS

Length: 77" // Width: 8.25" // Weight: 18.25 lbs.

CRASH CERTIFICATION

Fully meets SAE J3027

Ambulance Litter Integrity, Retention, and Patient Restraint

(Front: 30 mph, 22.5 G // Side: 15 mph 26G // Rear: 10mph/11G)

PRODUCT WARRANTY

3 Year Warranty



WHO IS Q'STRAIT?

Q'STRAIT began with a simple question:
How can wheelchair passengers be safely secured in vehicles?

Nearly 40 years ago, we answered with **the world's first fully integrated 4-point wheelchair passenger securement system**. Today, Q'STRAIT continues to build on this tradition of innovation, advancing its mission of providing safe accessibility to everyone by developing the automotive industry's highest quality wheelchair passenger safety solutions. As the globally recognized industry leader in wheelchair securements, we design, engineer and test our safety solutions for integration with nearly every major automotive OEM in the world.

Our mission is to **make safety accessible** for everyone.

And now, we are bringing that same innovation, engineering, passion and rigorous testing to emergency and non-emergency vehicles. Just like we believe that mobility passengers deserve access to the same level of safety as those around them we also believe that **First Responders, Caregivers and those who require Medical Transportation** should be as safe as possible in their most vulnerable moments.

Q'STRAIT®

4031 NE 12th Terrace
Oakland Park, FL 33334

Tel: 800-987-9987
Fax: 954-986-0021
Email: cs@qstraint.com



[QSTRAINT.COM/AMBULOK](https://qstraint.com/ambulok)

© 2023 Q'STRAIT®. All Rights Reserved. Q'STRAIT®, Q'STRAIT® logo, and AMBULOK® are registered trademarks of Q'STRAIT®, Inc.

MKM